


# Crono Super-PID Pump

For primary immunodeficiency

Based on the Crono standard pump, the **Super PID** pump incorporates more programming and safety features to provide facilities specifically for the subcutaneous infusion of immunoglobulins. They have a higher drive pressure setting to ensure trouble free infusion at comparatively high rates of flow.

The **Crono Super-PID** model can use dedicated 10 ml or 20 ml syringe reservoirs, and all are supplied with a range of carrying accessories for user convenience. The single battery cell has a typical life expectancy sufficient for 60 to 80 full syringe infusions.

Additional models are available in the *Crono S-PID* range which offer additional syringe capacities of 30, 50 and 100mL for greater volume administration. They also offer the choice of delivery program expressed in ml/hr or hours of infusion.

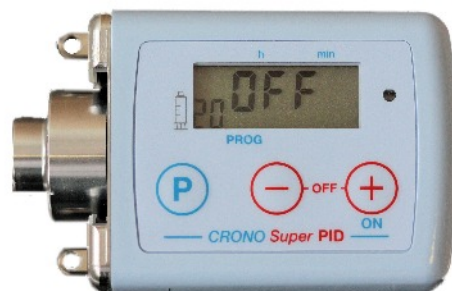

Specifications	
Dimensions/Weight	76x49x29mm/115g
Syringes (Dedicated)	10 or 20ml user selectable
Delivery times	1 to 99 hours in 15 min steps 15 to 60 mins in 5 min steps
Delivery Rate	N/A
Bolus Size	Fixed 0.2ml
Flow accuracy	+/- 2%
Keyboard Lock	2 security levels
Data Storage	All settings retained in memory even if battery removed
Accessories	
	Each pump is supplied with a range of accessories for patient convenience, which include: <ul style="list-style-type: none"> <li>• Soft Case</li> <li>• Battery cover removal tool</li> <li>• Soft waist belt</li> <li>• Neck Strap</li> <li>• Spare battery</li> </ul>

LR003/19BM

